

Highlights

Hot key for fast access
up to 10 saved programs

Unique electronic
tip ejection

Soft keys for
programming


Adjustment wheel

- Offers extremely fast volume setting and easy menu navigation
- Enables ergonomic one-handed operation
- Controls accurately manual pipetting and titrating speed with just a light touch of thumb

Operating button

Microwell plate tracker

- Unique built-in tracker for 96 and 384 well plates guides the user to pipette into correct wells
- Improves work efficiency and reliability of results
- Tracker is an optional feature used in conjunction with Pipetting, Reverse Pipetting and Multi Dispensing modes

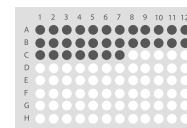

New generation technology

- Enhanced DC-motor concept takes Picus a step beyond the existing technological solutions providing outstanding performance
- Electronic brake stops the piston movement rapidly and precisely ensuring high precision especially in serial dispensing
- Optical sensor controls and monitors piston movement in real time ensuring unbeatable accuracy and reliability

